

PRODUCT FITTING


Gun Type X Ray Machine



Charging Base




Shielding Hood(Optional)

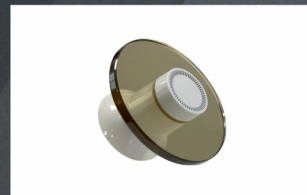

Power Supply




Plug

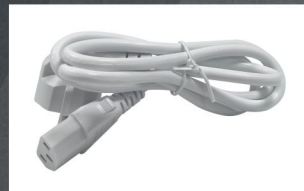

Remote Controller



User Manual

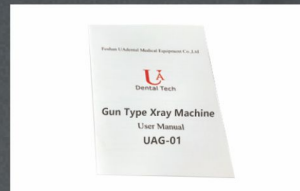

NEW LAUNCHED 2024

UAG-01

GUN TYPE X-RAY



OPTIONAL SHIELDING HOOD TO MEET
DIFFERENT REQUIREMENT

Foshan UAdental Medical Equipment Co.,Ltd

📍 Add: Building 2, South China International Medical Equipment Industry Center,
Shishan Town, Nanhai District, Foshan City, Guangdong Province, China



📞 Mob: +8613702756227

✉ Email: info@uadentech.com

📘 FB: Foshan UAdental

📷 Ins: foshanuadental

🌐 Web: www.uadentaltech.com

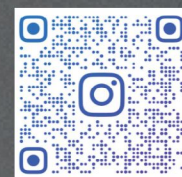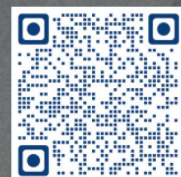

Foshan UAdental Medical Equipment Co.,Ltd

WHY UAG-01?


1. INDUCTION CHARGING BASE

The **first** Gun Type Xray with Induction Charging Base, Improve User Experience.




2. MILLISECOND EXPOSURE

65KV Tube Voltage, 2.5mA Tube Current for Millisecond Exposure , Effectively Avoid The Impact Of Trembling On Image Quality.

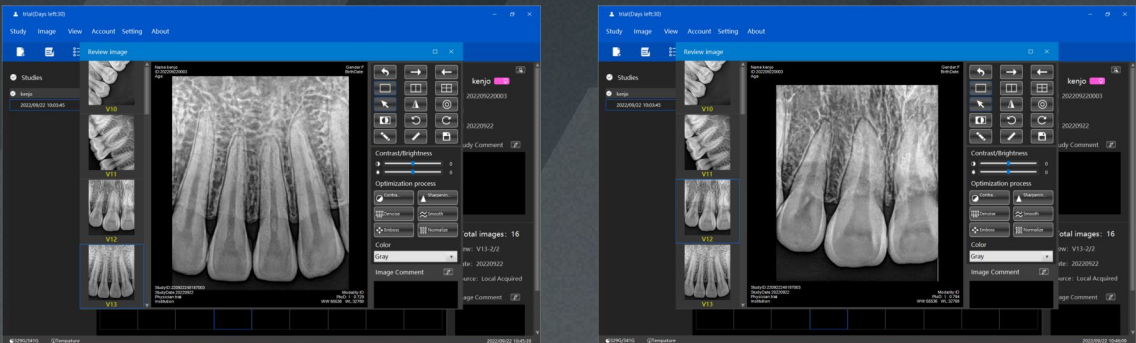

3. 0.4MM FOCUS TO GET HIGH DEFINITION IMAGES

High-voltage Generator With 0.4mm Focus, Enjoy Excellence Image Quality.

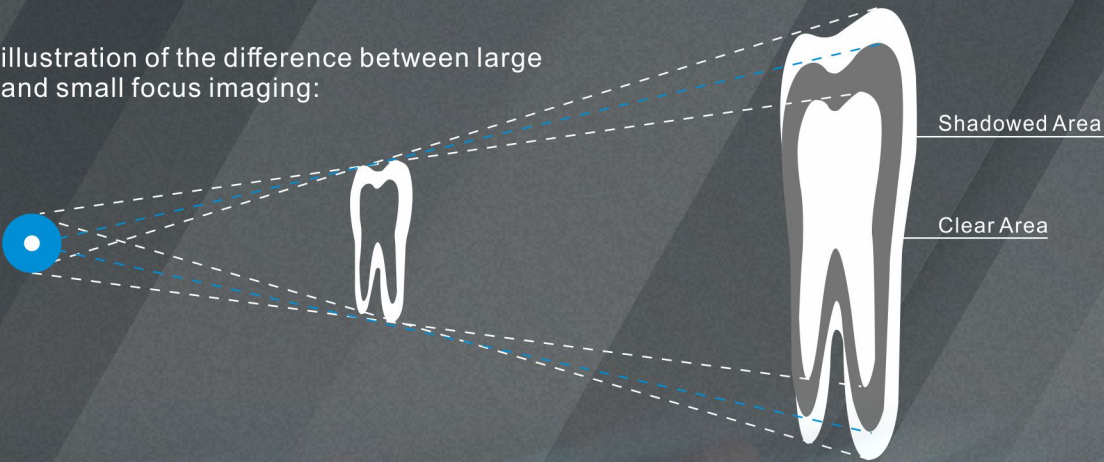

4. BIG CAPACITY RECHARGEABLE BATTERY—2500MAH

2500mAh X 4S High-capacity Lithium Battery, More Than 500 Shots Can Be Taken Continuously When The Battery Is Fully Charged.



5. INTELLIGENT ANGULATION SYSTEM

Real-Time Display Of Holding Angle, Shooting Angle Can Be Adjust As Patients' Posture

6. LOW RADIATION LEAKAGE AND LOCKING SYSTEM TO KEEP YOU SAFE

Leakage Radiation $\leq 0.05\text{mGy/h}$, Meet International Safety Standard, Fast Screen Lock To Prevent False Touch.

7. EASY OPERATION



Seven tooth position options



Three device options:
Sensor/Scanner/film



Two target groups:
Adult/Child



Press Both Buttons Together
To Open/Close The Child Lock



Press Both Buttons Together
To Enter System Setting

8. 2 YEARS REPLACEMENT WARRANTY

TEACHNICAL DATA

Model	UAG-01	Total Filtration	1.0mmAl
Tube Voltage	65kV	Battery Capacity	2500mAh
Tube Current	2.5mA	Maximum Inverter Frequency	80kHz
Power Supply	100-240VAC/50-60Hz	Weight	1.8kg
Exposure Time	0.01s-2s	Radiation Leakage Rate	1 Meter 0.05mGy/h
Continuous Shooting Times	>500 times	Warranty	2 years
Focus	0.4mm	Packing Size	40 x 37 x 23cm